



Domus Auto

Alternating Pressure Redistribution System With Automatic Adjustment


- Suitable for patients with medium to high risk of pressure ulcers
- Automatic adjustment for simple set up
- Cell-on-Cell bladder to avoid bottoming out
- Visual and audible alarm enhances patients safety and security






Domus Auto

Alternating Pressure
Redistribution
System With Automatic
Adjustment

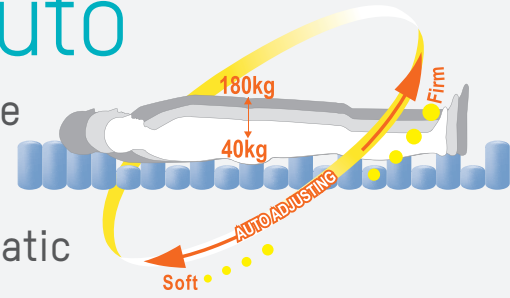

Cell-on-Cell Bladder
The bottom cells acts as basic pressure support to avoid bottoming out in case of power failure

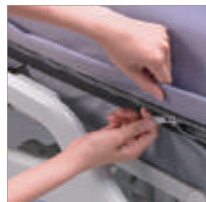

Zipper-all-Around Cover
Reduce ingress of unwanted fluid and simply one step to remove for cleaning



Automatic Adjustment
For simple set up



Tuning
To fine tune the pressure & optimize patient comfort



CPR Valve
Quickly deflate the mattress during emergencies



Nylon/PU Cover
Vapor permeable and water resistant with basic flame retardant level



Biocompatibility



Vapour Permeable



Water Resistance



Flame Resistance



Quick Connector
Once with cap on, instantly switch to static air mattress for transportation need




Visual and Audible Alarm

- Low Pressure
- Power Failure
- Tech Support


Pressure Relief During Alternating

Pressure mapping reveals that in alternating mode, bladder A & B take turns to inflate and deflate can significantly relieve the body pressure

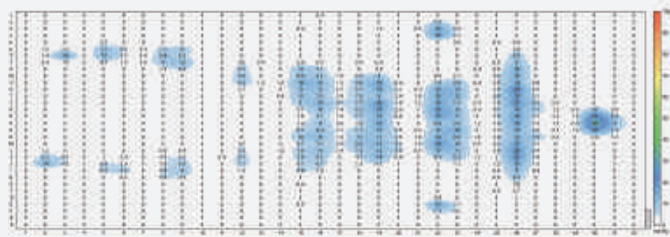

5 Modes to Fulfill Individual Needs

- Max Firm
- Alternating
- Continuous Low Pressure (CLP)
- Seat Inflation w/ Alternating
- Seat Inflation w/ CLP

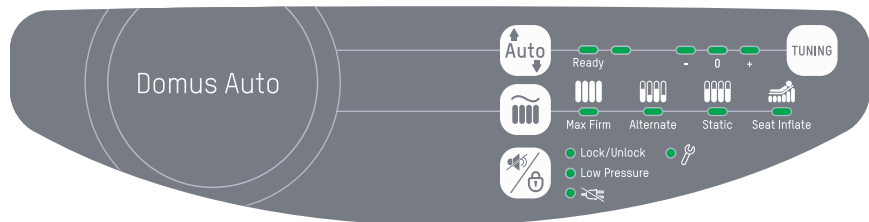

Specification

Model	Domus Auto	
Pump	Dimension	11.4" x 7.3" x 5.0" / 29 x 18.5 x 12.6 cm
	Weight	5.1 lb / 2.3 Kg
	Case Material	Flame retardant ABS
	Supply Voltage	AC 100-120V 60 Hz, 0.17A
	Operating Cycle	10 mins
Mattress		8" Replacement
	Dimension	78.7" x 35.4" x 8" / 200 x 90 x 20.3 cm
	Cell Height	8" x 20 cells
	Weight	15.2 lb / 6.9 Kg
	Material	Cover: Nylon/PU; Cell: Nylon/PU; Base: Nylon/PVC
	Maximum Patient Weight	450 lb / 200 Kg
	Flame Retardant Standards	CAL 117
HCPCS code	E0277	

Specifications are subject to change for improvement without notice

Ordering Information (For North America Only)

Model	Domus Auto	
System	With Quilted Cover	PD03501
	With Non-Quilted Cover	PD03502
Mattress	With Quilted Cover	PM09401
	With Non-Quilted Cover	PM09402
Pump	PD03001	

HEADQUARTERS

Wellell Inc.

No. 9, Minsheng St., Tucheng Dist.,
New Taipei City 236044, Taiwan
T +886 2 2268 5568
F +886 2 2268 6525
marketing@wellell.com

PAN-ASIA

CHINA

Wellell (Kunshan)Co., Ltd.
No. 1368, ZiZhu Road, Kunshan,
Jiang Su Province, China
T +86 512 8617 1660 Ext 1150
F +86 512 8617 8663
customer.cn@wellell.com

THAILAND

Apex Medical (Thailand)
No.111/152 Soi Pho Kao Yak 19,
Nawamin, Beung Kum,
Bangkok 10240,
Thailand.
T +66 02 018 2140
sales.thai@wellell.com

PAN-EUROPE

SPAIN

Apex Medical S.L.
Elcano 9, 6ª planta 48008
Bilbao. Vizcaya. Spain
T +34 94 4706408
F +34 94 4706409
info@apexmedical.es

UK

Wellell UK Limited
Unit 33, Great Western Business Park
McKenzie Way, Worcester ,WR4 9GN
United Kingdom
T +44 (0) 190 5774 695
sales.uk@wellell.com

FRANCE

Wellell France S.A.S
4 Boulevard de la Chanterie
49124 Saint-Barthelemy d' Anjou
France
T +33 2 72 79 27 13
F +33 2 85 73 07 67
info.france@wellell.com

GERMANY

SLK Vertriebsgesellschaft mbH
Am Herdicksbach 18, 45731 Waltrop
Germany
T +49 (0) 231-925360-0
F +49 (0) 231-925360-29
info@slk-gmbh.de
www.slk-gmbh.de

AMERICAS

USA

Wellell America Corp.
927 Mariner Street,
Brea, CA 92821, USA
T +1 714 671 3818
F +1 714 494 8173
sales.usa@wellell.com

Distributed by